

**GaleMed**

# Respiratory Care

Medical Technology For  
Your Health





# CONTENTS

7

Oxygen Nasal Cannula/  
Oxygen Mask

18

High Flow Nasal System

21

Bubble Humidifier & Large Volume  
Nebulizer

24

Aerosol Delivery/ Mask

34

GiO Manometer

39

Breathing Trainer

48

Product Index

51

Symbols and Cleaning Instructions





# Respiratory Care

Oxygen Nasal Cannula	08
Oxygen Tubing with Flow Indicator	08
Oxygen Conserving Nasal Cannula	09
Softip-1 Nasal Cannula	09
Flare Nasal Cannula	09
Oxygen Mask w/ Headgear	11
Simple Oxygen Mask	11
Face Tent	12
Oxi.M1 High-Oxygen Mask	12
OM-1/ OM-2/ OM-3 High-Oxygen Mask	13
Virus-Free Oxygen Mask	13
Tracheostomy Mask	14
Tracheal Tee Adaptor	14
M3 Mask	15
Capno Face Mask	16
Lite Adjustable Venturi Mask (LAV)	16
Fixed Venturi Diluters and Masks	17
High Flow Nasal Cannula System	19
Oxygen Tubing	20
Swivel Oxygen Tubing Connector	20
Disposable Flow Meter Adaptor	20
Bubble Humidifier	22
Large Volume Nebulizer	23
Aerosol Face Mask	25
Aerosol Face Tent	25
Tracheal Mask	25

## GaleMed

Oxygen Therapy is necessary for patients with low oxygen saturation who are unable to obtain adequate oxygen on their own. It consists of various interfaces for different requirements, such as low-flow oxygen devices, high-flow oxygen devices, and aerosol therapy.

GaleMed offers comprehensive solutions for patients with asthma, COPD, bronchopulmonary diseases, and dyspnea, integrating our four core values—safety, ergonomics, comfort, and practicality—into our product design.

### Making Every Breath Easier





Oxi.M1 Neb.Easy Mask	26
Neb.Easy Nebulizer/Neb.Easy Nebulizer Kit	27
Neb.Plus RF Nebulizer/ Neb.Plus RF Nebulizer Kit	28
Neb.Super High Efficiency Nebulizer	29
Neb-Valve Tee/	30
Regular Flow Nebulizer with Neb-Valve Tee Kit	
Neb.Plus LF1 Nebulizer/ Neb.Plus LF1 Nebulizer Kit	31
Dofin™ Spacer	32
MDI Valved Chamber Spacer	32
MDI Adaptor	33
Durable Spacer Mask/ Economic Spacer Mask	33
GiO™ Digital Pressure Gauge	35
GiO™ MIP / MEP Kit	36
GiO™ NIF Filtration Valve	37
Eco-Cuff Syringe Kit	37
Gas Sampling Line	38
ECHO Percussor	40
ECHO Expiratory Muscle Trainer	41
Tri-Ball Incentive Spirometer	42
Bravo/Dofin™ Breating Trainer	43
Expiratory Muscle Trainer	44
Peak Expiratory Flow Meter	44
Manual Percussor	45
Smile-Vac™ Manual Suction Unit	45
Mouth piece	46
Nose Clip	46
Filter for Spirometry	46



## Oxygen Nasal Cannula





## Oxygen Nasal Cannula

OXI.Plus™

## Oxygen Nasal Cannula

**Soft Ear Padding**

- Prevents pressure sores around ears
- 100% cotton material
- Minimizes skin irritation

**Magic Fitting**

- Expandable and flexible
- Easy for cleaning
- No adhesive needed

**Allergy-Free Silicone Prongs**

- Minimizes skin irritation
- Cured for better stability and flow path
- Ultra soft and lightweight
- Available in 4 sizes

- AO0040** Oxi.Plus™ Silicone Nasal Cannula Set, #0  
**AO0041** Oxi.Plus™ Silicone Nasal Cannula Set, #1  
**AO0042** Oxi.Plus™ Silicone Nasal Cannula Set, #2  
**AO0043** Oxi.Plus™ Silicone Nasal Cannula Set, #3

Package: 20 pcs/box

- Ultra soft and lightweight nasal cannula for oxygen therapy.
- Allergy-free silicone material to minimize skin irritation.
- Flexible and soft tube padding reduces pressure sores around the ears.
- Use with flow indicator to indicate and control airflow for patient safety.
- Expandable and detachable tubing for easy cleaning.
- Unique contoured nasal prongs to protect the philtrum and nares.
- Each nasal cannula loop includes 6mm ID oxygen tubing connector and Oxi.Plus nasal cannula (oxygen tubing with flow indicator available separately).
- For single patient use only.
- Conforms to ISO-10993 bio-compatibility standard.

**Oxygen Tubing with Flow Indicator**

- Indicates possible occlusions.
- Can detect airflow as low as 1 LPM.
- For single patient use only.

- AO0029** Flow Indicator w/Tubing, 160 cm  
**AO0030** Flow Indicator w/Tubing, 200 cm

Package: 25 pcs/box

OXI.Plus™

## Oxygen Nasal cannula


## Oxygen Nasal Cannula

**Oxygen Conserving Nasal Cannula**

- Patented dual lumen design enables accurate monitoring and sensing of patient's breathing pattern.
- Compatible with most oxygen conserving devices.

**AO0141** Oxygen Conserving Cannula

Package: 80 pcs/carton

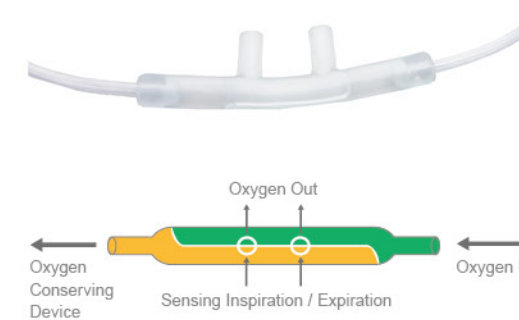**Standard Softip-1 Nasal Cannula**

- Includes nasal cannula and oxygen tubing (2-meter).
- Conforms to ISO-10993 biocompatibility standard.
- For single patient use only.

**Curved Softip-3 Nasal Cannula****Flare Nasal Cannula For Adult**

REF	Description	Type	QTY
<b>AO0066</b>	Softip-1 Nasal Cannula, Infant	Infant, Standard	100 pcs/carton
<b>AO0067</b>	Softip-1 Nasal Cannula, Child	Child, Standard	100 pcs/carton
<b>AO0068</b>	Softip-1 Nasal Cannula, Adult	Adult, Standard	100 pcs/carton
<b>AO0070</b>	Softip-3 Nasal Cannula, Infant	Infant, Curved	100 pcs/carton
<b>AO0071</b>	Softip-3 Nasal Cannula, Child	Child, Curved	100 pcs/carton
<b>AO0144</b>	Softip-3 Nasal Cannula, Adult	Adult, Curved	100 pcs/carton
<b>AO0072</b>	Softip-3 Nasal Cannula, Adult Large	Adult Large, Curved	100 pcs/carton
<b>AO0004</b>	Flare Nasal Cannula, Adult	Adult, Flare	100 pcs/carton







## Oxygen Mask



oxi.Plus™

Oxygen Mask


### Oxygen Mask w/ Headgear



- Phthalate-free, non PVC eco-friendly material.
- Ultra soft nose seal with sponge fitting reduces leaks and pressure sores.
- Two-way fixture that can be worn over the head or ears for greater comfort and stability.
- Sturdy mask shell resists distortion during transportation and storage.
- 55% lighter than conventional PVC O<sub>2</sub> masks.
- Each set includes a 2-meter oxygen tubing.
- 2-in-1 connector for attachment to oxygen tubing or nebulizer.
- For single patient use only.

**AO0115** Oxi.M1 Oxygen Mask Set, Size M


Package: 25 pcs/box

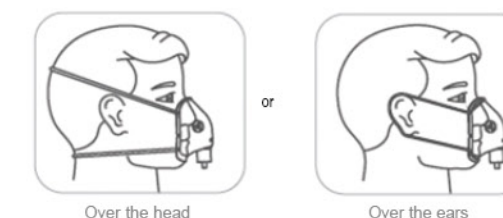

Over the head

Over the ears

### Simple Oxygen Mask

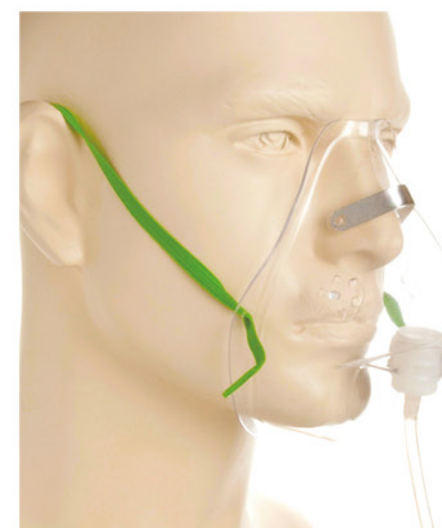

OM-2



OM-3

- Soft and lightweight mask designed for oxygen therapy.
- Aluminum nose clip for mask shape retention to reduce air leak around nose bridge.
- Clear and anatomic design for patient comfort and visual assessment.
- Latex-free elastic head strap for better fit and less skin irritation.
- Each set includes 2-meter oxygen tubing and oxygen tubing connector.
- Conforms to ISO-10993 biocompatibility standard.
- For single patient use only.

**AO0086** OM-2 Oxygen Mask Set, Adult

**AO0085** OM-2 Oxygen Mask Set, Child

**AO0084** OM-2 Oxygen Mask, Adult

**AO0005** OM-2 Oxygen Mask, Child

**AO0153** OM-3 Oxygen Mask Set, Child

**AO0154** OM-3 Oxygen Mask Set, Adult

Package: 50 pcs/box



## Face Tent



- Designed for patients with facial skin burn or trauma.
- Conforms to ISO-10993 biocompatibility standard.

**AO0112** OM-5 Face Tent with Oxygen Tubing

**AO0111** OM-5 Face Tent with Oxygen Adaptor (no tubing)

**Package:** 25 pcs/box


## Oxi.M1 High-Oxygen Mask



- Provides oxygen concentration levels between 85% - 99%.
- One-way valve between the mask and reservoir ensures reservoir oxygen concentration and minimizes CO<sub>2</sub> retention.
- Phthalate-free, non PVC eco-friendly material.
- Ultra soft nose seal with sponge fitting reduces leaks and pressure sores.
- Two-way fixture that can be worn over the head or ears for greater comfort and stability.
- Sturdy mask shell resists distortion during transportation and storage.
- 55% lighter than conventional PVC O<sub>2</sub> masks.
- For single patient use only.

**AO0039** Oxi.M1 Hi-Oxygen Mask Set, Size L

**Package:** 50 pcs/box

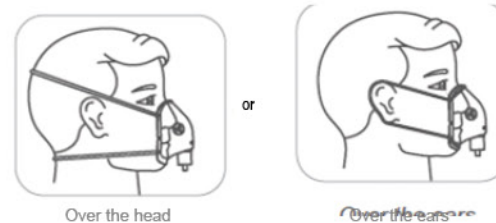

## High-Oxygen Mask



OM-1



OM-2

OM-3

- Provides oxygen concentration levels between 85% - 99%.
- One-way valve between the mask and reservoir ensures reservoir oxygen concentration and minimizes CO<sub>2</sub> retention.
- Ergonomic mask shape fits a wide range of face sizes.
- OM-2/OM-3: Flexible and malleable aluminum nose clip ensures a secure fit and reduces air leaks around the nose bridge.
- Each set includes a Hi-Oxygen Mask, and a 2-meter oxygen tubing with regular oxygen tubing connector.
- Conforms to ISO-10993 biocompatibility standard.
- For single patient use only.

**AO0089** OM-1 Hi-Oxygen Mask Set, Adult

**AO0088** OM-1 Hi-Oxygen Mask Set, Child

**AO0093** OM-2 Hi-Oxygen Mask Set, Adult

**AO0094** OM-2 Hi-Oxygen Mask Set, Child

**AO0095** OM-3 Hi-Oxygen Mask Set, Adult

**AO0097** OM-3 Hi-Oxygen Mask Set, Child

**AO0096** OM-3 HiO<sub>2</sub> Mask, Adult w/Tubing, PVC Bag

**Package:** 50 sets/box

## Virus-Free Oxygen Mask



- Designed for infection control of contagious airborne diseases such as SARS and Bird-Flu.
- Includes Airsoft™ mask and silicone headgear with superior seal to minimize oxygen leak and increase patient comfort.
- Filters bacteria and viruses up to 99.999%.
- Features an oxygen reservoir for high oxygen concentration.
- One-way valve reduces CO<sub>2</sub> retention.
- Swivel elbow design for flexible use.
- Made with medical grade, latex-free material to minimize allergies and other complications.
- For single patient use only.

**AO0164** Virus-Free Oxygen Mask Child

**Package:** 25 sets/box



## Tracheostomy Mask



- Contoured mask design to provide better seal.
- Flexible 360° double-swivel elbow connector with suction port.
- Includes a mask and a 2-meter oxygen tubing with regular plus connector.
- Conforms to ISO-10993 biocompatibility standard.

**AO0106** OM-6 Trachea Mask, Adult w/Tubing  
**AO0105** OM-4 Trachea Mask, Adult w/Tubing  
**AO0104** OM-4 Trachea Mask, Child w/Tubing

**Package:** 50 sets/box

## Tracheal Tee Adaptor



- Oxygen reservoir features 15 cm tube for temporary oxygen storage.
- Oxygen supply port with cap: tapered 6mm.
- 22/15 mm and 22M ISO standard connection ports.
- Each set includes 22M, 22/15 mm Tracheal Tee piece and 15cm disposable corrugated tubing.
- Conforms to ISO-10993 biocompatibility standard.

**AO0058** Trachea Tee Set, w/Tubing

**Package:** 100 sets/box

## M3 Mask



The M3 masks give oxygen therapy a whole new meaning, providing safety and comfort where it is needed. Instead of using PVC (Poly-vinyl-chloride), M3 masks opt for a safer and more environmentally friendly material to create well-fitted, hard-shelled masks with comfortable cushion that fit different facial features.

### M3 Mask

- Fitting cushion gently seals around the nose bridge, preventing gas flow from damaging the eyes.
- Reduces the back pressure during exhalation.
- Rigid Non PVC and clear shell resists deformation during shipment.

**AO0024** M3 Oxygen Masks, L  
**AO0025** M3 Oxygen Masks, M

**Package:** 25 set/box

The M3 masks give oxygen therapy a whole new meaning, providing more flexibility with 360° swivel oxygen tubing connector.




#### Sturdy Oxygen Tubing (M3 Oxygen Masks only)

- Reinforced lumen for crush-resistant tubing.




#### Soft and Elastic Head Strap

- Ultra-comfortable material.
- Can be worn around the ears or over the head.

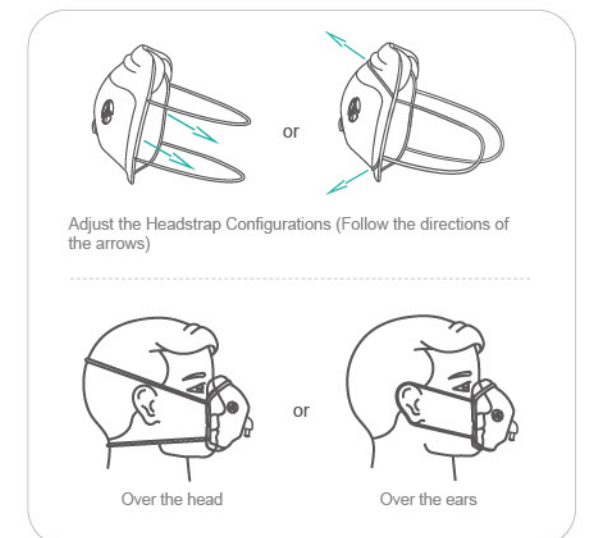



## Capno Face Mask

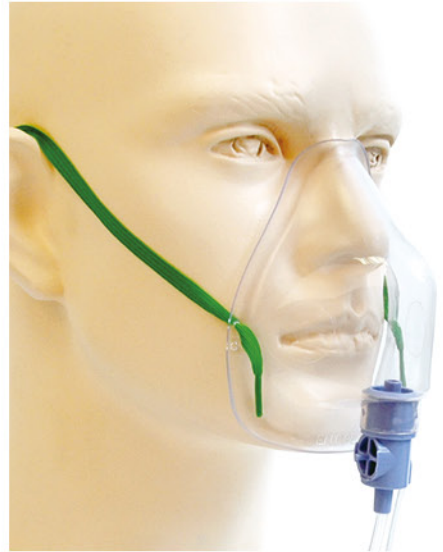

- Designed for capnography.
- Supports both nasal and oral gas exchange measurements.
- Compact design that combines CO<sub>2</sub> sampling tubing, oxygen tubing and exhalation outlet.
- Anatomic mask design provides better seal for accurate CO<sub>2</sub> sampling.
- Features flexible 360° swivel connector.
- 4.3 Luer Lock for end-tidal CO<sub>2</sub> device connection.
- 4.8 mm ID O<sub>2</sub> connector.
- Includes 2-meter oxygen tubing.

**AO0003** Capno-Mask, Adult

**AO0002** Capno-Mask, Child

Package: 10 sets/box

## Lite Adjustable Venturi Mask (LAV)



- Designed based on Bernoulli's Principle.
- 11 different oxygen concentration settings. Provides oxygen concentration levels between 24% - 90% (at 3 LPM).
- Compact and lightweight.
- Features a convenient pull-twist-lock mechanism to reduce the risk of accidental mis-adjustments.
- Venturi adaptor with 22M connector for mask, face tent and other applications.
- Each set includes mask, oxygen tubing and adjustable oxygen diluter.
- Can be used with trachea mask, trachea tee and face tent.
- For single-patient use only.
- Conforms to ISO-10993 biocompatibility standard.

**AO0045** LAV Venturi Mask Set, Adult

**AO0162** LAV Venturi Mask Set, Child

**AO0163** LAV Venturi Mask Set, Infant

Package: 10 sets/box

**AO0001** Adjustable Oxygen Diluter

**AO0044** Lite Adjustable Venturi Adaptor, 22F

Package: 100 pcs/carton



AO0044



AO0001

## Fixed Venturi Diluters and Masks

- Designed based on Bernoulli's Principle.
- Provides oxygen concentration levels between 24% - 60%.
- Available in 7 different settings: 24, 28, 31, 35, 40, 50 and 60%.
- Color coded for easy identification.
- 22M or 22F connection port.
- For single patient use only.
- Conforms to ISO-10993 biocompatibility standard.



• With 22M Connector

SAP	Color	LPM	Oxygen
<b>AO0118</b>	Blue	4	24%
<b>AO0119</b>	Yellow	4	28%
<b>AO0120</b>	White	4	31%
<b>AO0056</b>	Green	6	35%
<b>AO0057</b>	Red	8	40%
<b>AO0121</b>	Orange	12	50%

Package: 100 pkgs/box

• With 22F Connector

SAP	Color	LPM	Oxygen
<b>AO0048</b>	Blue	2	24%
<b>AO0049</b>	Yellow	4	28%
<b>AO0050</b>	White	4	31%
<b>AO0051</b>	Green	6	35%
<b>AO0052</b>	Red	8	40%
<b>AO0053</b>	Orange	12	50%

Package: 100 pkgs/box



### Fixed Venturi Face Mask

- Each set includes 6 fixed venturi diluters with an oxygen face mask, a venturi hood and a 2-meter oxygen tubing.
- 22M Adaptor for connection to nebulizer.

**AO0098** AM-1 Venturi Mask Set, Adult

**AO0099** AM-3 Venturi Mask Set, Adult

**AO0102** AM-3 Venturi Mask Set, Adult

**AO0101** AM-2 Venturi Mask Set, Child

Package: 50 pcs/box



### Fixed Venturi Tracheostomy Mask

- Each set includes 6 fixed venturi diluters with a tracheal mask, a venturi hood and a 2-meter oxygen tubing.
- 22M Adaptor for connection to nebulizer.
- Swivel elbow connector for greater flexibility.

**AO0100** AM-4 Venturi Tracheal Mask Set, Adult

Package: 50 pcs/box





# High Flow Nasal System



OXI.Plus™

High Flow Nasal System

## High Flow Nasal Cannula System




### NHF Prong Kit

#### Advanced oxygen management

The High Flow Nasal Cannula is designed to deliver high flow via ultra soft silicone nasal cannula to improve patient breathing conditions and comfort. This High Flow Nasal Cannula is available to combine with our Humi. AIDE humidification technology which offers an effective and efficient solution for oxygen therapy.

- AO0009 NHF Prong Kit, #2
- AO0010 NHF Prong Kit, #3
- AO0011 NHF Prong Kit, #4

Package: 100pcs/carton

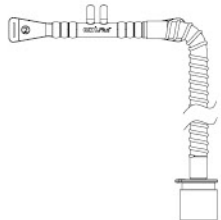

### High Flow Nasal Cannula



#### A Softer Touch

High Flow Nasal Cannula is made of super soft silicone material and encompasses a special undercut design which mimics the anatomical curve of philtrum (the area under nose). The anatomical curve design reduces the friction and irritation on skin and helps stabilize the cannula on the patient as well as safely and comfortably accommodate patient nares, allowing for optimal oxygen delivery.

#### Quick Release


The "snap-on" button on one end of the cannula is designed to enable easy removal of the headgear or adjustment.

#### Condensation Drainage

Accumulated in the nasal cannula can be quickly drained via one of the two tubing connection ports, without any interruption to the air supply.

#### Sizes Available

High Flow Nasal Cannula is a patient interface for high flow oxygen therapy. It features a simple yet intuitive contoured design, lending itself naturally to improve comfort and functional usability. High Flow Nasal Cannula comes in 3 sizes, catering to a variety of patients with different sizes of nasal passages.

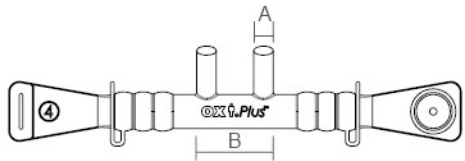

Spec.	Size 2	Size 3	Size 4
Diameter ( A )	4.8	5.5	6.0
Distance ( B )	4.2	6.5	10

Unit: mm



# Oxygen tubing & Connector

oxi.Plus™

## Oxygen Tubing



- Kink-free tubing features an inner star lumen design to ensure continuous oxygen flow.
- Two styles of oxygen tubing connector available (C22 Regular Plus and C23 Fits-all).
- Fits-all connector for flow meter, humidifier or other oxygen/air supply devices.
- Different lengths and types of oxygen connectors available upon request.

REF	Description	QTY
AO0062	C23 Star Lumen Tubing, 200 cm, fits-all connector	100 pcs/box
AO0063	C22 Star Lumen Tubing, 200 cm, 6F	100 pcs/box
AO0064	C22 Star Lumen Tubing, 180 cm, 6F	100 pcs/box
AO0065	C22 Star Lumen Tubing, 1000 cm, 6F	30 pcs/box

## Swivel Oxygen Tubing Connector



- Tapered 6M connection on both ends.
  - Also available with clamp.
- AO0046 Swivel Connector for Oxygen Tube  
AO0047 Swivel Connector w/Clamp

Package: 100 pcs/box

## Disposable Flow Meter Adaptor



- Available with a variety of connection types at flow meter outlets.
- AO0074 Flow Meter Connector - Extended (Green)  
AO0075 Flow Meter Connector (Green)  
AO0076 Flow Meter Connector (White)

Package: 50 pcs/pkg, 500 pcs/box




## Bubble Humidifier & Large Volume Nebulizer

...



## Bubble Humidifier ⓘ

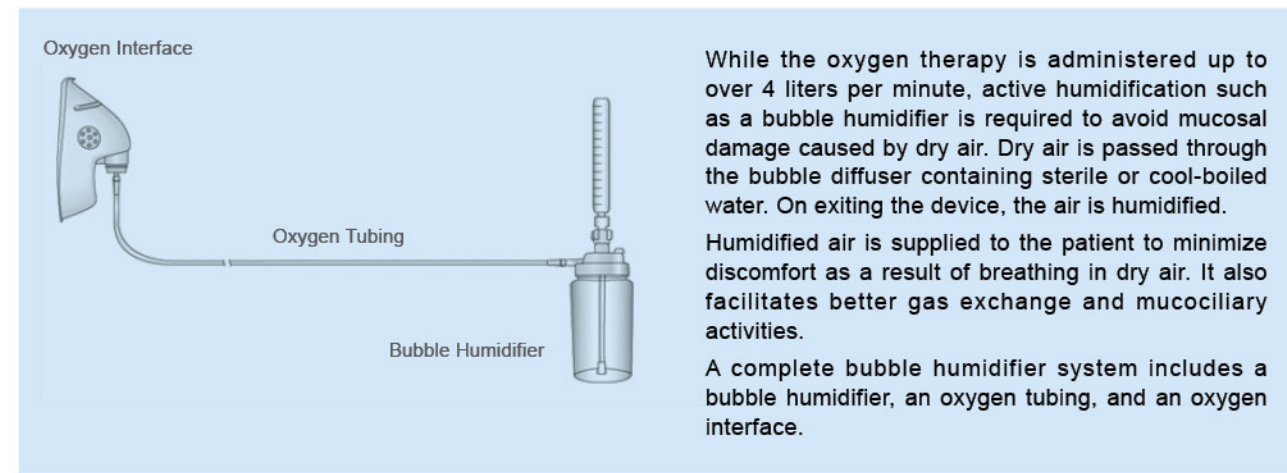

## BH-1 Bubble Humidifier



## Features:

- We provide a range of pressure safety valves which relieve the pressure at the levels of 2, 4, 6 PSI.
- Volume: 400 ml.

## BH-2 Bubble Humidifier



## Features:

- Consists of specially-designed connecting hole for heater rod to avoid unnecessary risk of overheating.
- Safety Valve: 4 PSI.
- Volume: 400 ml.

## BH-3 Bubble Humidifier




## Features:

- Smaller size of bubble humidifier for oxygen-supply devices and for short-term usage.
- Safety Valve: 4 PSI.
- Volume: 180 ml.

REF	Description	Pressure Safety Valve	QTY
AH0028	BH-1 Bubble Humidifier, 2 PSI	2 PSI	24 pcs/carton
AH0030	BH-1 Bubble Humidifier, 4 PSI	4 PSI	24 pcs/carton
AH0035	BH-1 Bubble Humidifier, 6 PSI	6 PSI	24 pcs/carton

REF	Description	Type	QTY
AH0027	BH-3 Bubble Humidifier	Smaller Volume	24 pcs/carton
AH0034	BH-2 Bubble Humidifier	Heater Rod Insertion	24 pcs/carton

## Large Volume Nebulizer ⓘ

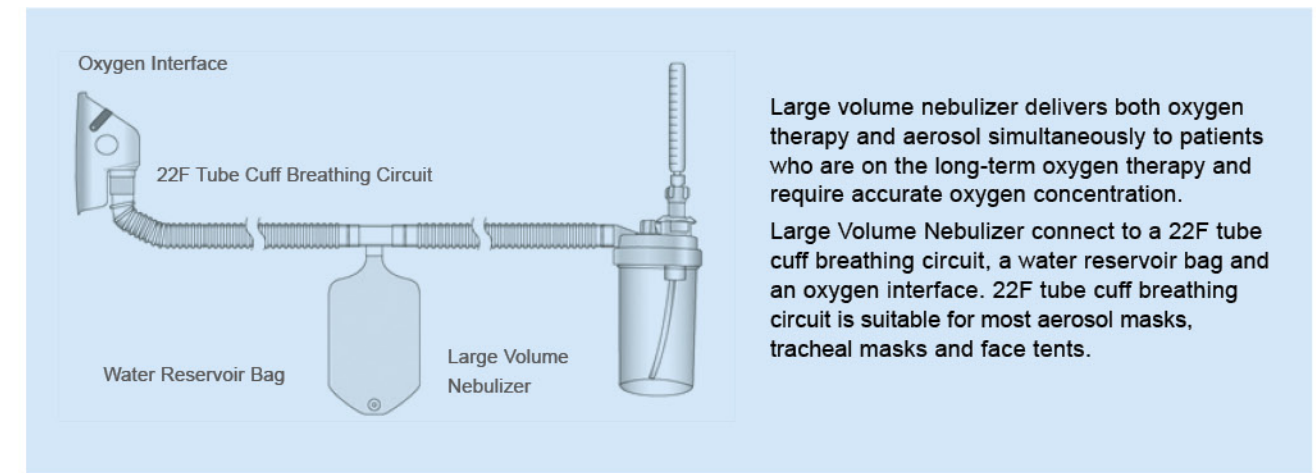

## Large Volume Nebulizer



## Features:

- Venturi-Style Entrainment Ring: Adjust the oxygen concentration with incremental dial according to patient needs.
- Heater Rod Insertion Port: Allow heater rod insertion to deliver warm humidified aerosol.
- Oxygen concentration range from 35% to 100%.
- Set flow at 10 l/min for output of 60 ml/hr.
- Volume: 400 ml.

REF	Description	QTY
AH0032	LVN Nebulizer, PP	24 pcs/carton
AB0101	Water Reservoir, 600 ml	50 pcs/carton
AT0012	EVA Hose, 22F, 180 cm	50 pcs/carton
AT0001	EVA Hose, 22F, 100FT	1 pc/carton

## Water Reservoir Bag



## 22F EVA Breathing Circuit







# Aerosol Delivery



## Large Volume Nebulizer System Interface

### Aerosol Face Mask



Shown with adult type for demonstration only

**AD0064** AM-1 Aerosol Mask, Infant

Package: 50 pcs/carton



**AD0068** AM-2 Aerosol Mask, Child

**AD0069** AM-2 Aerosol Mask, Adult

Package: 50 pcs/carton



**AD0067** AM-3 Aerosol Mask, Adult

Package: 50 pcs/carton

### Aerosol Face Tent



**AO0110** AM-5 Face Tent

Package: 25 pcs/ctn

### Tracheal Mask



**AO0108** AM-6 Trachea Mask, Adult

**AO0103** AM-4 Trachea Mask, Child

Package: 50 pcs/ctn



## Neb.Easy Nebulizer




### Oxi.M1 Neb.Easy Mask

- Phthalate-free, non PVC eco-friendly material.
- Ultra soft nose seal with sponge fitting reduces leaks and pressure sores.
- Two-way fixture that can be worn over the head or ears for greater comfort and stability.
- Lower residual drug in the reservoir and delivers fine particles.
- Each set includes Oxi.M1 aerosol mask and Neb. Easy Nebulizer Set.

**AD0026** Oxi.M1 Neb.Easy Mask Set, Size M


**AD0027** Oxi.M1 Neb.Easy Mask Set, Size L

Package: 25 pcs/ctn

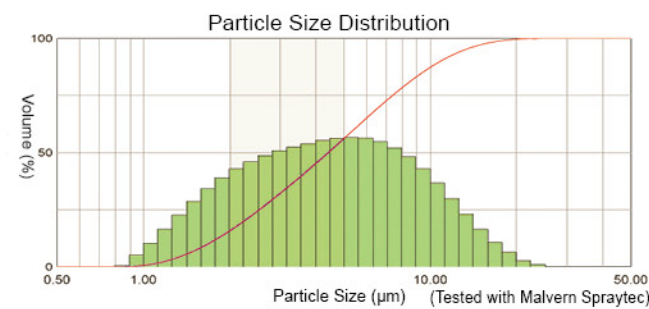

Particle Size(MMAD)	Gas Flow Rate	Reservoir Volume
3.8 µm	0.25 ml/min	6 ml

## Neb.Easy Nebulizer

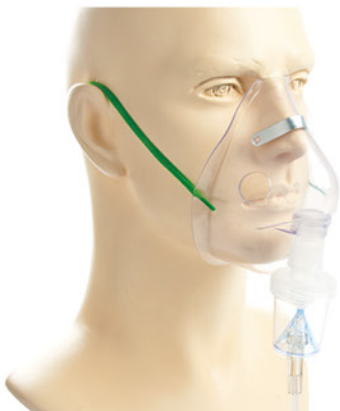

**AD0054** Neb.Easy AM-2 Mask, Adult w/Tubing

**AD0058** Neb.Easy AM-2 Mask, Child w/Tubing

Package: 50 pcs/ctn



**AD0002** Neb.Easy Nebulizer, 18M

**AD0096** Neb.Easy Nebulizer, 22F

Package: 50 pcs/ctn

## Neb.Easy Nebulizer Kit



- With 45° L-shape Mouth Piece.
- Patients could hold the device in lower position, which makes the inhalation therapy more comfortable.

**AD0042** Neb.Easy Nebulizer Kit & 45° L-shape Mouth Piece

Package: 50 pcs/box

## Nebulizer Kit



- Neb.Easy Nebulizer with 18M or 22F connector.

**AD0044** Nebulizer Kit (Nebulizer with 18M Connector)

**AD0040** Nebulizer Kit (Nebulizer with 22F Connector)

Package: 50 pcs/ctn

## Isolated Nebulizer Kit



- Tapered 6M connection on both ends.
- Also available with clamp.

**AD0035** Isolated Nebulizer Kit with standard cuff oxygen tubing

**AD0036** Isolated Nebulizer Kit with fits-all cuff oxygen tubing

Package: 25 pcs/ctn



## Neb.Plus RF Nebulizer




### Neb.Plus Regular Flow Nebulizer

- Equipped with a port for medication re-fill.
- Can be connected to a breathing circuit to enable drug delivery during mechanical ventilation.
- V-shape Reservoir optimizes medication and minimizes medical residue.
- Suitable for aerosolized medication at gas flows of 6 to 10 LPM.

**AD0071** Neb-1 Nebulizer Mask, Adult

**AD0070** Neb-1 Nebulizer Mask, Child

**Package:** 50 pcs/ctn

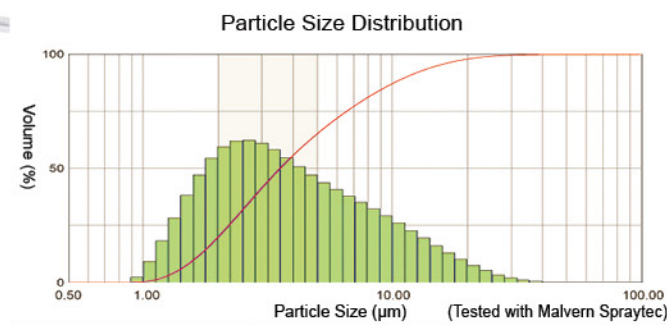

Particle Size(MMAD)	Gas Flow Rate	Reservoir Volume
3.5 µm	0.3 ml/min	20 ml



**AD0037** Neb-1 RF1 Nebulizer Kit

**Package:** 25 pcs/ctn



**AD0022** RF1 Nebulizer Kit with T Adaptor 22F, 22/15 mm

**Package:** 25 pcs/ctn

## Neb.Super High Efficiency Nebulizer




### Neb.Super Mask Kit

- Reduces average treatment time by up to 45%.
- Effective aerosol delivery system allows for enhanced gas exchange through medical deposition of the surrounding respiratory tract and lungs.
- Suitable for aerosol therapy at gas flows of 6 to 10 LPM.
- Each set includes aerosol mask and Neb.Super Nebulizer Set.

**AD0019** Neb.Super Mask Kit, Child

**AD0020** Neb.Super Mask Kit, Adult

**Package:** 25 pcs/ctn

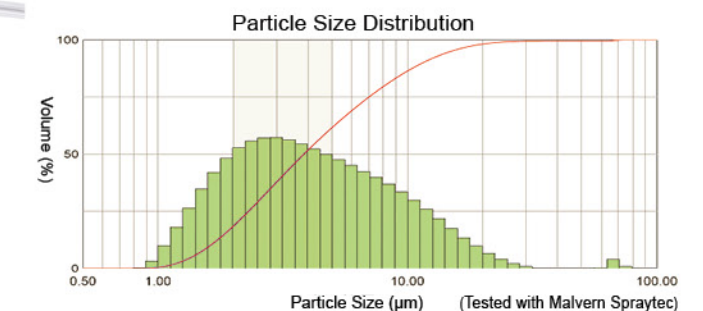

Particle Size(MMAD)	Gas Flow Rate	Reservoir Volume
3.7 µm	0.63 ml/min	20 ml



**AD0018** Neb-Super Nebulizer Kit

**Package:** 25 pcs/ctn




## Aerosol Delivery



## Neb-Valve Tee

The Neb-Valve Tee employs a T-shaped valve for use in support of a nebulizer within a breathing circuit. With the use of the Neb-Valve Tee, patients are able to receive aerosol therapy without interrupting ventilation. The Neb-Valve Tee can be used with or without a nebulizer; the valve automatically closes itself when a nebulizer is not used and provides only ventilated gas. However, when a nebulizer is locked into the valve, the mechanism allows the valve to open to deliver both gas and aerosol medication.

- Pressure & Flow Rate: < 0.29 cmH<sub>2</sub>O @ 60 LPM (AV0076, AV0097, AV0077)
- < 0.50 cmH<sub>2</sub>O @ 60 LPM (AV0075)
- < 3.65 cmH<sub>2</sub>O 60 LPM (AV0074)

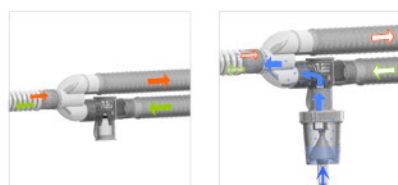

AV0076 / AV0097 / AV0077



AV0075



AV0074

REF	Description	Pressure & Flow Rate	QTY
AV0076	Valve Tee, 22F, 18F, 22/15mm	< 0.29 cmH <sub>2</sub> O @ 60 LPM	50 pcs/box, 4 box/ctn
AV0097	Valve Tee, 22F, 18F, 22M	< 0.29 cmH <sub>2</sub> O @ 60 LPM	50 pcs/box, 4 box/ctn
AV0077	Valve Tee, 22M, 18F 22M	< 0.29 cmH <sub>2</sub> O @ 60 LPM	50 pcs/box, 4 box/ctn
AV0075	Valve Tee, 15F, 18F, 15M	< 0.50 cmH <sub>2</sub> O @ 60 LPM	50 pcs/box, 4 box/ctn
AV0074	Valve Tee, 10F, 18F, 10M	< 3.65 cmH <sub>2</sub> O 60 LPM	50 pcs/box, 4 box/ctn

## Regular Flow Nebulizer with Neb-Valve Tee Kit



- Regular flow nebulizer at 6-10 LPM for use with adult ventilators.
- Luer medicine filling port can be directly attached to drug delivery pump or syringe.
- Each set consists of RF1 nebulizer, Neb-Valve Tee and 160 m gas supply tubing.

- AD0004 RF1 Nebulizer w/ATV, 1.6M Tube
- AD0005 RF1 Nebulizer w/Valve T Kit, 22F, 22M
- AD0006 RF1 Nebulizer w/Valve T Kit, 22M, 22M

Package: 25 pcs/ctn

## Neb.Plus Low Flow Nebulizer



AD0016 Neb.Plus LF1 Nebulizer

Package: 25 pcs/box



- Low flow nebulizer at 1-3 LPM can reduce the impact on neonatal ventilators and flow measurement.
- Luer medicine filling port can be directly attached to drug delivery pump or syringe.
- Each set consists of LF1 nebulizer, Neb-Valve Tee and 160m gas supply tubing.

AD0023 Neb.Plus LF1 Nebulizer Kit with Neb.T Valve 22F, 22/15

Package: 25pcs/ ctn



- Nebulizer flow rate: 1-3.5 LPM.
- MMAD: 2.03µm/1 LPM.
- Nebulization rate: ≈ 0.02 ml/min.
- Nebulizer Volume: 20 ml.
- Medication fill port for long term therapy.
- DEHP-Free and Latex-Free material to minimize the risk of skin irritation.

AD0017 Neb.Plus LF2 Nebulizer w/ Supplemental Gas Line

Package: 25 pcs/box



## Aerosol Delivery

### Dofin™ Spacer

- Enhances drug delivery rate by up to 25% and reduces medical waste.
- Adaptor compatible with most existing MDI in the market.
- Optimization of the medicated air travel dynamics.
- Lightweight and easy to grasp and hold.



Metered Dose Inhalers (MDI) are used to treat respiratory problems through aerosol therapy. One of the major disadvantages of MDI's is the need for coordination between MDI actuation and patient inhalation. In most cases, patients fail to use MDI's properly. The use of a spacer can greatly help patients breathe at their own pace and reduce harmful side effects caused by inadequate drug deposition or trapping in the oropharynx and large conducting airways. Using an MDI Spacer allows the user to inhale medication straight into the lungs and optimizes treatment efficiency.




**AD0079** Dofin™ Spacer w/Mouth Piece

Package: 48 box/carton



**AD0080** Dofin™ Spacer for Breathing Circuit

Package: 18 sets/carton

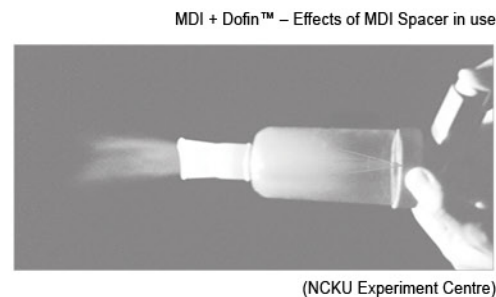

- AD0074** MDI Valved Chamber Spacer, 22M  
**AD0076** Spacer with Silicone Mask #0  
**AD0077** Spacer with Silicone Mask #1  
**AD0078** Spacer with Silicone Mask #2

Package: 48 pcs/ carton

## MDI Adaptor



- Allows metered dose inhaler (MDI) medication directly deliver to the ventilator circuit.
- Universal port accommodates both plastic and metal MDI canisters. Available with adult and pediatric circuit sizes.
- MDI canister stem is inserted directly into port of MDI adaptor.
- Cap preserves closed circuit integrity.

**AD0031** MDI Adaptor, 22F, 22B

**AD0032** MDI Adaptor, 15F, 15B

Package: 50 pcs/carton

**AD0033** MDI L-Adaptor, 22M/15F, 15M

Package: 100 pcs/carton

## Spacer Mask



- AM0089** Durable Silicone Spacer Mask #2  
**AM0090** Durable Silicone Spacer Mask #3  
**AM0091** Durable Silicone Spacer Mask #4


Package: 48 pcs/carton



- AM0087** Economic PVC Spacer Mask (Child)  
**AM0088** Economic PVC Spacer Mask (Adult)

Package: 100 pcs/carton



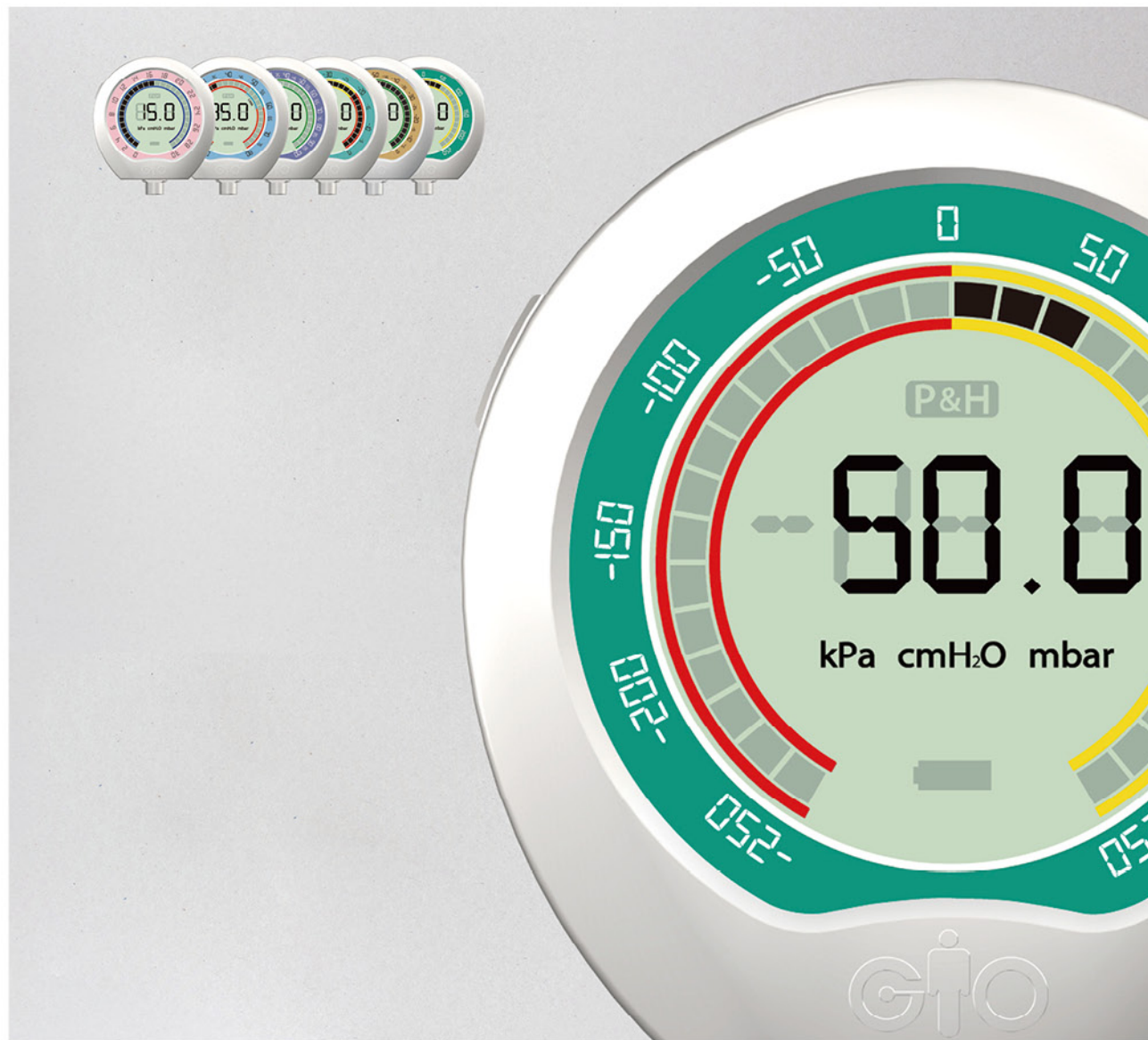

## GiO™ Manometer



## GiO™ Manometer

### GiO™ Digital Pressure Gauge



GiO™ Digital Pressure Gauge is an economical, portable and compact device enables the measurement and display of pressure values in a single instrument to fulfill pressure monitoring needs. Besides, we have different specifications of gauge to meet the needs of different situations.



- 1 GiO™ Manometer (AQ0015 / GiO™ Manometer, -250~250cmH<sub>2</sub>O)
- 2 Elbow Connector 7.6M, 6B
- 3 Straight Luer 7.6M, 4.3M Leur
- 4 Elbow Connector 7.6M, 4.3B
- 5 Hydrophobic Filter 7.6M, 4.3F Leur
- 6 Pressure line 50cm, 4.3M Leur
- 7 Silicone Pressure Line 50cm, 4.3M Leur
- 8 Lithium Battery

REF	Description	QTY
AQ0008	GiO Manometer, 0~30 cmH <sub>2</sub> O	10 pcs/carton
AQ0012	GiO Manometer, 0~80 cmH <sub>2</sub> O	10 pcs/carton
AQ0013	GiO Manometer, -30~100 cmH <sub>2</sub> O	10 pcs/carton
AQ0014	GiO Manometer, -60~0 cmH <sub>2</sub> O	10 pcs/carton
AQ0015	GiO Manometer, -250~250 cmH <sub>2</sub> O	10 pcs/carton
AQ0016	GiO Manometer, -100~0 cmH <sub>2</sub> O	10 pcs/carton
AF0003	Hydrophobic Filter, 7.6M w/Luer	10 pcs/carton




## GiO™ Digital Pressure Gauge Accessories

Where  
GiO  
Pressure  
Gauge  
is used



### Maximal Inspiratory and Expiratory Pressure (MIP) & (MEP)

The MIP and MEP are pressure measurement values that indicate the respiratory muscle strength of the user. Pressure gauge records the peak inspiratory and expiratory pressure. Understanding patient MIP and MEP is important as it can help evaluate voluntary breathing capability, determine ventilatory capacity and can also indicate the possible development of respiratory symptoms. Furthermore, MIP & MEP measurement is also helpful in evaluating the success of weaning patients from mechanical ventilators. For healthy individuals, physical activities result in an increase in oxygen demand where the higher the intensity of exercise, the more oxygen is needed. Respiratory muscle strength is highly related to oxygen intake and hence respiratory muscle training is very important. MIP or MEP assessment is a convenient method to evaluate lung function and provides clear assessment of training results.

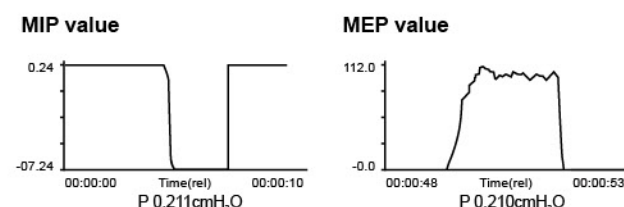

References:

1. ATS/ERS Statement on Respiratory Muscle Testing.
2. Rocha Crispino Santos, M.A. Pinto, M.L. Sant'Anna, C. & Bernhoeft, M. 2010 Maximal respiratory pressures among adolescent swimmers. For more references, please visit our website at [www.gio-solutions.com](http://www.gio-solutions.com)
- 2010 Meodimal respiratory pressures among adolescent swimmers



### GiO™ MIP / MEP Kit

- Connect to **GiO™**(AP0005) for MIP or MEP measurement.
- Each kit includes 2 separate One Way Valve (MEP valve and MIP valve), bacterial filter and mouthpiece.
- A soft, anatomic and comfortable mouthpiece provides better seal.
- Single use interface minimizes the risk of cross contamination.

- 1 GiO™ Manometer (AQ0015 / GiO™ Manometer, -250~250cmH<sub>2</sub>O)
- 2 Hydrophobic Filter
- 3 Silicone Pressure Line
- 4 One Way Valve
- 5 Bacterial Filter with 4.3F Port
- 6 Mouthpiece



REF	Description	QTY
AP0005	MIP/MEP Kit	20 pcs/carton

## GiO™ Digital Pressure Gauge Accessories

### Negative Inspiratory Force (NIF)


Ventilator Dependent Patient

Calculating the respiratory muscle strength of mechanically ventilated patients who cannot sustain spontaneous breathing can be complicated. Measuring chronic ventilator-dependent patients' NIF is the first step towards weaning from the ventilator.

#### The NIF Filtration Valve

The NIF Filtration Valve is an innovative interface used with GiO™ to provide NIF readings. From the concept of a flute, user can occlude it to prevent air from escaping during inhalation. The built-in filter prevents the moisture and secretion from damaging the pressure gauge and provides the optimal measurement condition.

### Why measure NIF

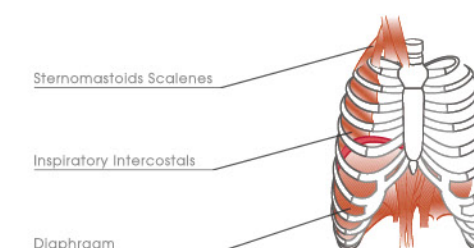

#### Muscles of respiration

NIF indicates the strength of inspiratory muscles, helping clinicians to evaluate whether a patient could be weaned from mechanical ventilation. The value of NIF also provides reliable reference of the improvement.



- 1 GiO™ Manometer (AQ0014 / GiO™ Manometer, -60~0 cmH<sub>2</sub>O)
- 2 NIF Valve With Filter
- 3 Collapsible Tube
- 4 Mouthpiece

### GiO™ NIF Filtration Valve

- Single use interface minimizes the risk of cross contamination.
- Flexible extension tube with standard 15mm connector.
- Filtration valve regulates and filters air and sputum during measurement.
- Mouthpiece or connectors for tracheostomy tube and endotracheal attachment.

AP0003 NIF Kit

Package: 30 pcs/carton



- 1 GiO™ Manometer (AQ0012 / GiO™ Manometer, 0~80 cmH<sub>2</sub>O)
- 2 Hydrophobic Filter
- 3 3-way stopcock
- 4 Pressure line
- 5 Luer lock Syringe

### Eco-Cuff Syringe Kit

- Eco-Cuff is designed to inflate and monitor the cuff pressure of endotracheal tube or tracheal tube.
- Single use interface minimizes risk of cross-contamination.

AP0001 Eco-Cuff Syringe Kit

Package: 5 pcs/carton



## Gas Sampling Line



Different types of tubes for delivering or monitoring airflow or CO<sub>2</sub>.

- AB0005** Gas Sampling Line, 300 cm, 4M
- AB0060** Gas Sampling Line, 300 cm, 4F, 4M
- AB0061** Gas Sampling Line, 50 cm, 4M, 4M
- AB0062** Gas Sampling Line, 300 cm, 4F, 4F
- AB0064** Gas Sampling Line, 180 cm, 4M, 4M

Package: 50 pcs/carton



## Breathing Trainer




## ECHO Percussor

ECHO Percussor is a unique device which could not only clean sputum in the respiratory tract but also improve the expiratory muscle strength and cough efficacy by using positive airway pressure and vibration.



- Provides desired airway oscillation frequency to mobilize mucus and strengthen expiratory muscle.
- Unique quick adjustment knob to vary oscillation frequency or level for lung training.
- Indicator for airflow passing through the device to ensure proper use of device.
- Quick "Tit and Escape" design allows air be released or breathed in quickly to maximize its safety.

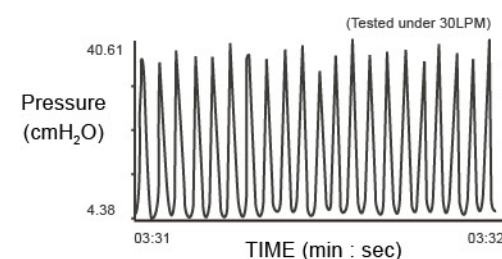

Adjustment Valve



- AW0061** ECHO Percussor
- AW0062** ECHO Percussor Kit, w/Tracheal Adaptor
- AW0063** ECHO Percussor w/ Tracheostomy Mask Size #0

Package: 60 sets/carton


Optimum Oscillation Frequency for Mucus Mobilization: 20~25 Hz

Using the ECHO Percussor is an easy and simple way to manage sputum retention in contrast to invasive methods such as suction. Echo Percussor helps loosen mucous surrounding the respiratory tracts, making airway clearance easier.




With Mouthpiece

- ISO 22mm Mouthpiece and Flex tube
- Can be used with a bacterial filter

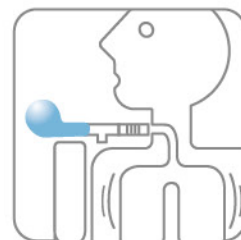

With Tracheostomy Tube


- ISO 22/15mm valved T Connector
- For use on tracheostomy patients
- Can be used with a bacterial filter

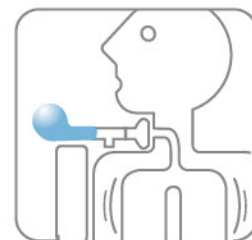

With Tracheostomy Mask

- ISO 22/15mm valved T Connector
- For use on tracheostomy patients
- Can be used with a bacterial filter

## ECHO Expiratory Muscle Trainer

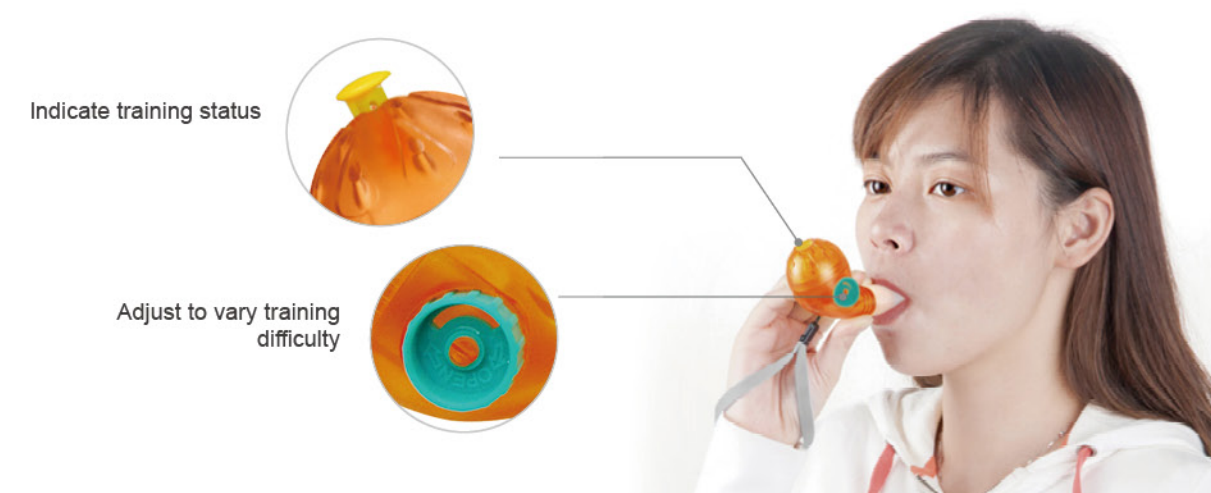

- An expiratory muscle trainer encourages users to blow through the device to lift the training indicator.
- Suitable for users to improve performances requiring maintenance and managing expiratory muscles.
- Adjustable pressure valve adjust training levels.
- Comes with a pop-up training status indicator when the weight has been lifted.
- Includes a safe-release neck strap and a soft, comfortable mouthpiece.



With Mouthpiece

- AL0009** ECHO Expiratory Muscle Trainer

Package: 24 pcs/carton

- ISO 22mm Mouthpiece and Flex tube
- Can be used with a bacterial filter



## Tri-Ball Incentive Spirometer



- Designed to stimulate lung function recovery after surgery.
- Breathing in through the spirometer helps strengthening the inspiratory muscles.
- Color-coded balls in each chamber to provide a visual incentive to indicate the inspiratory flow rate at 600 ml/sec, 900 ml/sec and 1200 ml/sec.

**AL0014** Tri-ball Incentive Spirometer


Package: 20 sets/carton

## Bravo/Dofin™ Breathing Trainer



- All-in-one trainer, for inspiratory and expiratory muscle training.
- More than 10 training levels for different resistance training.
- Each coloured trainer indicates different resistance level.
- Yellow ball designed to indicate breathing status and strength.
- Compact design, simple operation, portable, lightweight and easy to clean.

REF	Description	QTY
AL0019	Dofin™ Trainer, Blue	15 sets/carton
AL0027	Dofin™ Trainer, Blue	15 sets/carton
AL0020	Bravo™ Trainer, Green	32 sets/carton
AL0021	Bravo™ Trainer, Orange	32 sets/carton

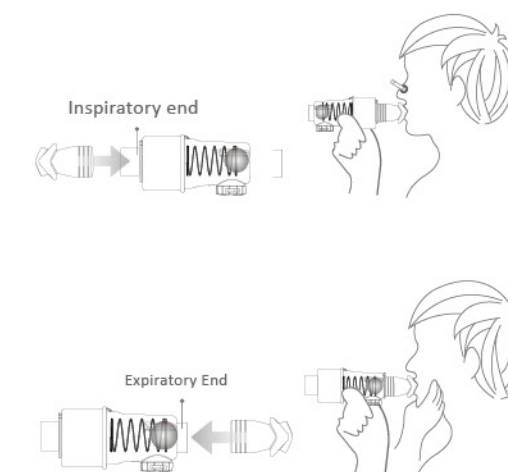

## Inspiratory/ Expiratory Resistance Values

Maximal Inspiratory Resistance Level			Maximal Expiratory Resistance Level		
Blue	Green	Orange	Blue	Green	Orange
39	79	163	33	82	154





## Expiratory Muscle Trainer



- Adjustable PEEP valve (5~20 cmH<sub>2</sub>O) allows patients to start with lower resistance then increase the level of resistance to strengthen their muscle.
- Optimize bronchodilation when combined with respiratory drug delivery via nebulizer.
- Improves strength, endurance and walking distance for patients with chronic obstructive pulmonary disease (COPD).
- Convenient to use. Easy to clean and assemble.

**AL0003** EMT Expiratory Muscle Trainer

Package: 24 boxes/carton

## Peak Flow Meter and Accessorie

- Measures the peak expiratory flow (PEF) score before and after medication for patients.
- Slip-free design with finger placement grooves.
- Lightweight & compact for travel.
- Visual indicators, slider, marker and record graphs are available for tracking PEF score easily.
- Simple self-assessment tables can review PEF score.
- Available in 2 models: 200~800 LPM for adults and 50~300 LPM for children.

**AL0012** Peak Flow Meter, 800 LPM  
**AL0013** Peak Flow Meter, 300 LPM

Package: 50 pcs/carton



Slip-Free Finger Grooves



Protective Seal



Hand Strap



## Manual Percussor



Manual percussor is an airway management device assisting chest physiotherapy. Manual percussor offers a convenient alternative to manual chest percussion via "hand cupping" with unique manual percussor cups, allowing caregiver to perform chest percussion with manual pounding or clapping to loosen secretions.

**AW0056** Manual Percussor, S  
**AW0057** Manual Percussor, M  
**AW0055** Manual Percussor, L

Package: 48 pcs/carton



1. Cup the manual percussor between your thumb and index finger as shown.
2. As patient leans forward, pat the manual percussor at a steady but firm rhythm.

## Smile-Vac™ Manual Suction Unit



- Manual suction device removes excessive mucus, saliva or liquid in the oral cavity or airway.
- Soft and bio-compatible suction tube minimizes tissue irritation.
- For single-patient use in homecare or in emergency.
- Easily clean.

**AW0014** Smile-Vac™ Manual Suction Unit

Package: 24 sets/carton



## Accessories for Spirometry and Lung Training

Package: 100 pcs/carton



## Nose Clip



**AL0010** Nose Clip, Blue  
**AL0011** Nose Clip, Orange

Package: 100 pcs/carton

## Filter for Spirometry




- Compatible with ISO 30 mm.
- Resistance: < 0.8 cmH<sub>2</sub>O @40 L/Min.
- Dead Space: < 100 ml.

**AL0022** 30 mm ID, 30 mm OD

Package: 100 pcs/carton

Users should fully understand the uses of products and the following symbols before use.

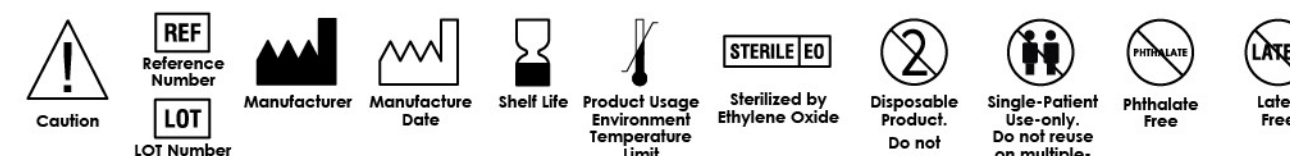

### Common Materials and Sterilization Method

Material	Autoclavable	Chemical Disinfection	Ethylene Oxide Sterilization
PVC	No	No	No
Acrylic	No	No	No
K-Resin	No	No	No
ABS	No	No	Yes
Polyethylene(PE)	No	Yes	No
Polypropylene(PP)	No	Yes	No
Hytrel	Yes (121c)	No	Yes
Polycarbonate(PC)	Yes (121c)	No	Yes
Polysulfone(PSF)	Yes	Yes	Yes
Silicone Rubber	Yes	Yes	Yes
Stainless Steel	Yes	Yes	Yes

**!** Disposable and Single-patient use products are designed to reduce the risk of cross-contamination and therefore have strict prohibitions on re-processing or reusing them on multiple users.

When re-using durable products, make sure to follow validated institution policy and procedure as well as the information available for disinfection processes.

#### When to clean or disinfect:

1. If product has been used
2. If product is contaminated
3. If product has been used for more than 24 hours on a single-patient

#### Manufacturer's suggested procedure for disinfection:

The following instruction is a general guide. Please refer to individual organization's operating instructions on cleaning and disinfecting products:

##### a. High Temperature Pressure Disinfection:

Products should be disinfected at 134°C for 5 minutes (eg. for materials such as polysulfone, silicone, stainless steel). For hytrel and polycarbonate products, the temperature should not exceed 121°C and should not be processed for longer than 15 minutes.

##### b. Chemical Disinfection:

Products should be soaked under warm water with mild detergents. Wash and rinse thoroughly to remove any residual material and store in a cool dry place to allow it to air dry.

##### c. Ethylene Oxide Disinfection:

Follow the instruction for use for Ethylene Oxide disinfection devices.

### Cleaning instructions for homecare products:

When cleaning or disinfecting homecare products, users can follow the following general guidance as well as the individual product user manuals:

1. Boiling Water can be used on products that are not degradable. Products that are metal and stainless steel nature can be disinfected by boiling for 10 - 15 minutes with water.
2. Hot-water disinfection (Pasteurization) is a high level, non toxic disinfection process that can be used to reprocess products. In this process, a previously cleaned equipment is immersed in a water bath at 63-65°C for 30 minutes.
3. Alcohol Disinfection: Place products under a combination of mild soap and water and rinse thoroughly, then; proceed to soak under 70% Isopropanol (Alcohol) for 5 minutes which then should be rinsed with sterile water.
4. Detergents: Clean and remove all organic matter (blood and tissue) along with other residue from product or equipment. Wash products using a mild detergent such as (glutaraldehyde) and then rinse with sterile water. Products should then be dried and stored in a cool and dry place.

#### Important Note:

1. After disinfection and sterilization procedures, products and equipment should be stored in an appropriate package so as to avoid further contamination
2. Disposable products should be disposed after one-time use. This is to minimize cross-contamination and ensure patient safety.



## Product Index

REF	Old REF	Page
AB0005	30002	38
AB0060	30001	38
AB0061	30005	38
AB0062	30009	38
AB0064	30016	38
AB0101	70030	23
AD0002	3551	26
AD0004	NV31	31
AD0005	NV32	31
AD0006	NV33	31
AD0016	NB01	30
AD0017	NB02	31
AD0018	NH51	29
AD0019	NM51	29
AD0020	NM52	29
AD0022	NT31	28
AD0023	NV11	31
AD0026	OM65	26
AD0027	OM66	26
AD0028	3001	46
AD0029	3004-A01	46
AD0030	3005	46
AD0031	30341	33
AD0032	30342	33
AD0033	30351	33
AD0035	3507	27
AD0036	3509	27
AD0037	3511	28
AD0040	3555	27
AD0042	3556	27
AD0044	3558	27
AD0054	3569	26
AD0058	3570	26
AD0064	3711	25
AD0067	3723	25
AD0068	3732	25
AD0069	3733	25
AD0070	3736	28
AD0071	3737	28
AD0074	4200	32
AD0076	4201	32
AD0077	4202	32
AD0078	4203	32
AD0079	4205	32

REF	Old REF	Page
AD0080	4206	32
AD0096	3553	26
AD0108	3739	28
AF0003	39015	35
AH0027	3402	22
AH0028	3403	22
AH0030	3404	22
AH0032	3406	23
AH0034	3407	22
AH0035	3416	22
AL0001	3002	46
AL0002	3004	46
AL0003	5481	44
AL0009	DK21	41
AL0010	3055	46
AL0011	3056	46
AL0012	4121	44
AL0013	4131	44
AL0014	4311	42
AL0019	DT11	43
AL0020	DT12	43
AL0021	DT13	43
AL0022	4410	46
AL0027	DT14	43
AM0087	5253	33
AM0088	5254	33
AM0089	5256	33
AM0090	5257	33
AM0091	5258	33
AO0001	33001	16
AO0002	3076	16
AO0003	3077	16
AO0004	3280	9
AO0005	3642	11
AO0009	NH22	19
AO0010	NH23	19
AO0011	NH24	19
AO0024	AME0007	15
AO0025	AME0008	15
AO0029	OF21	8
AO0030	OF22	8
AO0039	OM36	12
AO0040	ON10	8
AO0041	ON11	8

## Product Index

REF	Old REF	Page
AO0042	ON12	8
AO0043	ON13	8
AO0044	OV30	16
AO0045	OV31	16
AO0046	30019	20
AO0047	30022	20
AO0048	30120	17
AO0049	30121	17
AO0050	30122	17
AO0051	30123	17
AO0052	30124	17
AO0053	30125	17
AO0056	3084	17
AO0057	3085	17
AO0058	3096	14
AO0062	31220	20
AO0063	31310	20
AO0064	31510	20
AO0065	31610	20
AO0066	3231	9
AO0067	3232	9
AO0068	3233	9
AO0070	3240	9
AO0071	3241	9
AO0072	3243	9
AO0074	33114	20
AO0075	33115	20
AO0076	33116	20
AO0084	3643	11
AO0085	3644	11
AO0086	3645	11
AO0088	3657	13
AO0089	3658	13
AO0093	3675	13
AO0094	3678	13
AO0095	3679	13
AO0097	3680	13
AO0098	3683	17
AO0099	—	17
AO0100	3715	17
AO0101	3753	17
AO0102	3754	17
AO0103	5201	25
AO0104	5203	14

REF	Old REF	Page
AO0105	5204	14
AO0106	5206	14
AO0108	5208	25
AO0110	5231	25
AO0111	5239	12
AO0112	5240	12
AO0115	OM15	11
AO0118	3081	17
AO0119	3082	17
AO0120	3083	17
AO0121	3086	17
AO0141	3072	9
AO0144	3242	9
AO0153	3637-MIT	11
AO0154	3638-MIT	11
AO0162	OV32	16
AO0163	OV33	16
AO0164	3652	13
AP0001	GC22-01	37
AP0003	GN10	37
AP0005	GN30	36
AQ0008	GB10	35
AQ0012	GB20	35
AQ0013	GB30	35
AQ0014	GB40	35
AQ0015	GB60	35
AQ0016	GB50	35
AT0001	6001	23
AT0012	38002	23
AV0074	7VLD11	30
AV0075	7VLD21	30
AV0076	7VLD31	30
AV0077	7VLD33	30
AV0097	7VLD32	30
AW0014	1411	45
AW0055	39202	45
AW0056	39200	45
AW0057	39201	45
AW0061	DK11	40
AW0062	DK12	40
AW0063	DK13	40





Date     /     /

business card

Clean & Sterilization

MATERIALS USED FOR GALEMED PRODUCTS:


Single patient use or disposable items are prohibited to be reused.  
Reusable items are treated as below:

Material	Autoclave 134°C	EtO gas	Liquid Sterilizer
ABS	Distort	Yes	Yes
Polycarbonate (PC)	Distort	Yes	Excellent
Polyethylene (PE )	Marginal to poor	Yes	Excellent
Polypropylene (PP)	Yes	Yes	Excellent
Polysulfone	Excellent	Yes	Excellent
PVC	Marginal to poor	Yes	Distort
Silicone rubber	Yes	Yes	Excellent
Stainless, Aluminum	Excellent	Yes	Excellent
Rubber	Distort	Yes	Yes


CLEANING & STERILIZATION STEPS:

- 1 Disassemble the product as per the IFU, if necessary.
- 2 Wash components thoroughly in lukewarm water with a mild detergent which shall be compatible with the material of the device.
- 3 Rinse the components thoroughly in lukewarm water to remove all traces of detergent.
- 4 Sterilize the components using one of the following methods:  
a. Autoclaved (max 134°C /273°F).  
b. Ethylene oxide.  
c. The most common disinfecting solution is used the equipment coming from the contact with the patient. Rinse thoroughly with water after sterilization.
- 5 Dry all components thoroughly after sterilization.
- 6 Inspect all components prior to usage; if any tear or any damage, please replace it.
- 7 Repack properly after usage to avoid any contamination.


Symbol Guide For This Catalog




**Do not reuse**  
To indicate that the product is for single use only.




Indicates that natural rubber latex was not used in the manufacturing of the product, its container, or its packaging.



Sterilizable in a steam sterilizer (autoclave) at temperature 121°C (250°F) will be accepted



Indicates that product that does not contain the phthalate plasticizer DEHP.



Sterilizable in a steam sterilizer (autoclave) at temperature 134°C (273°F) will be sterilizable




**GaleMed** specializes in the design and manufacture of consumable medical devices in the field of respiratory care. Our products cater to the complete patient journey in respiratory care, including neonatal care, emergency care, anesthesia, critical care, and pulmonary rehabilitation. Backed by a team of clinical professionals, we provide customers with the breadth and depth of product solutions.

**GaleMed** is headquartered in Taiwan and operates manufacturing facilities in Taiwan and Xiamen. With distinct production lines and dedicated marketing teams in both locations, GaleMed has been committed to customers and market needs globally, spanning Europe, America, Southeast Asia, East Asia, and other parts of the world.

With more than 35 years of OEM and OBM experience, GaleMed provides one-stop services, including **Mold Design, Professional Process, Process Validation & Qualification, and Automatic Storage Retrieval System**. GaleMed possesses excellent quality control and global logistics capabilities.

**GaleMed** has an experienced RQA team and ensures that its employees remain properly certified and continuously trained. We possess medical device market approvals and quality system management certifications from multiple countries to keep pace with evolving medical device regulations, including **TFDA, NMPA, FDA, EU MDR, TW-QMS, ISO 13485, MDSAP, JGMP** and others. Thus, our products and services are in line with the latest international standards, helping maintain our customers in a competitive market.

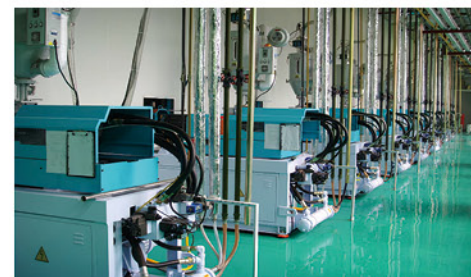

**Established in**  
1986

**Business Types**  
ODM/ OEM manufacturer,  
Contract Manufacturer & Exporter

#### Main Products

- Respiratory Care
  - Critical Care
  - Anesthesia Care
  - Emergency Care
  - Oxygen Therapy
  - Aerosol Therapy
  - Neonatal Care
  - Lung Training

#### Quality Management System

- ISO 13485
- ISO 9001
- MDSAP
- JAPAN J-GMP
- TAIWAN QMS

#### Product Certification

- EU MDR Certificate
- FDA Listing and 510K
- China CFDA Registration
- Taiwan TFDA Registration





**GaleMed**

**GaleMed**

**China GaleMed  
1<sup>st</sup> Facility  
(Xiamen)**

**China GaleMed  
2<sup>nd</sup> Facility  
(Xiamen)**

**China GaleMed  
Aplus Facility & Sales  
Department (Jimei)**

**Taiwan GaleMed  
Facility (Yilan)**

**GaleMed Corporation  
HQ (Taipei)**



**Research,  
Development  
and Manufacture**

Since 1986

**GaleMed Corporation**

[www.galemed.com](http://www.galemed.com)

Global Sales Office

**GaleMed Limited Taipei Branch**

5F., No. 209, Sec. 1, Fuxing S. Rd., Da'an Dist.,  
Taipei City 106, Taiwan (R.O.C.)

Email: [sales@galemed.com](mailto:sales@galemed.com)

Tel: +886 2 8791-9348



GaleMed



LinkedIn



FB